




Discontinuing an Infusor™ from a PICC Line

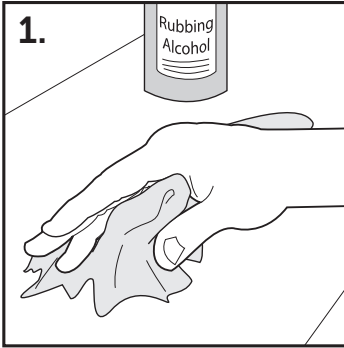

1. Clean working surface with rubbing alcohol.

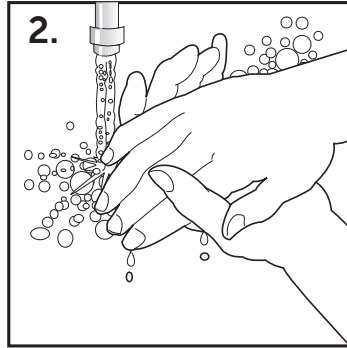


2. Wash hands and dry with clean towel.

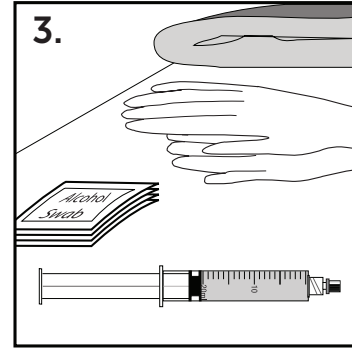


3. Gather equipment & place on clean surface:

- 20ml normal saline syringe
- Gloves
- Alcohol swabs
- Heparin*


*Your nurse will inform you if your PICC requires heparin.
*If your line **does** require heparin, repeat steps 7, 8, 9

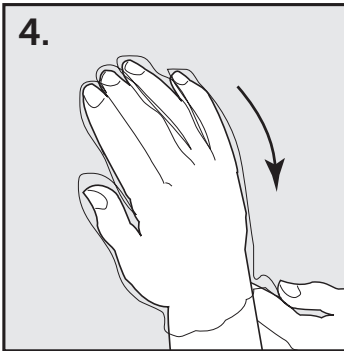


4. Put on gloves.

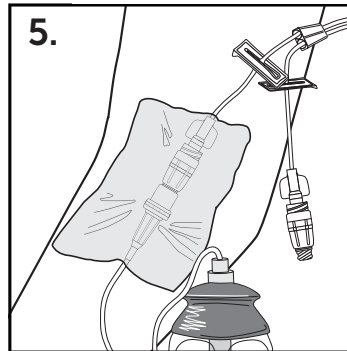


5. Slide clamp of line connected to infusor to the closed position.

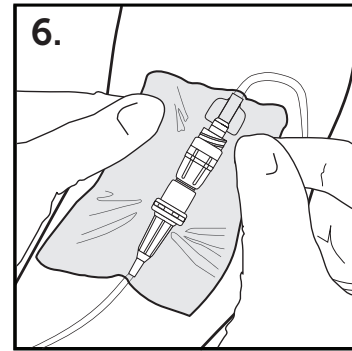


6. Remove the tape that holds white sensor against skin.

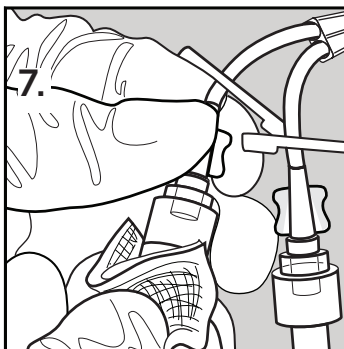

7. Scrub the cap of the line that is not connected to the infusor pump. Allow cap to dry.

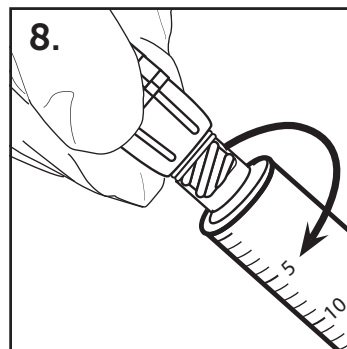


8. Connect the syringe of saline to the clean cap by twisting the syringe clockwise.

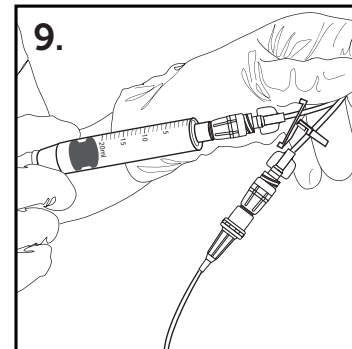


9. Open slide clamp and push 20mLs saline using a start-stop motion. While pushing, clamp the line when 1ml saline left.

□ If your PICC line requires heparin please repeat steps 7-9 with heparin syringe

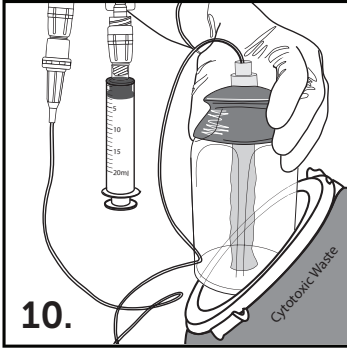


Disconnect tubing from PICC line and put in the Cytotoxic Waste container.

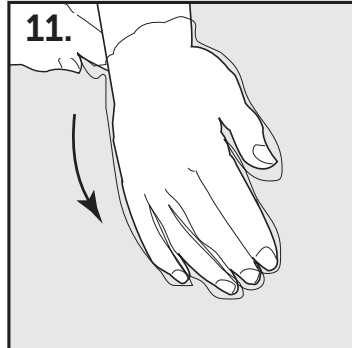

Take off the gloves and discard in cytotoxic waste container.



Clean working surface with rubbing alcohol.

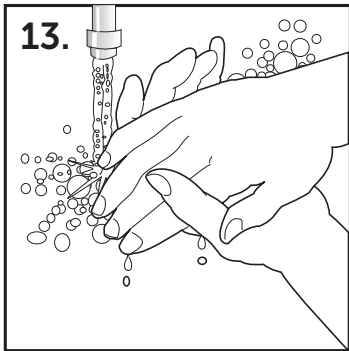

Wash hands and dry with clean towel.